


# JAW CRUSHER / BUCKET RSD

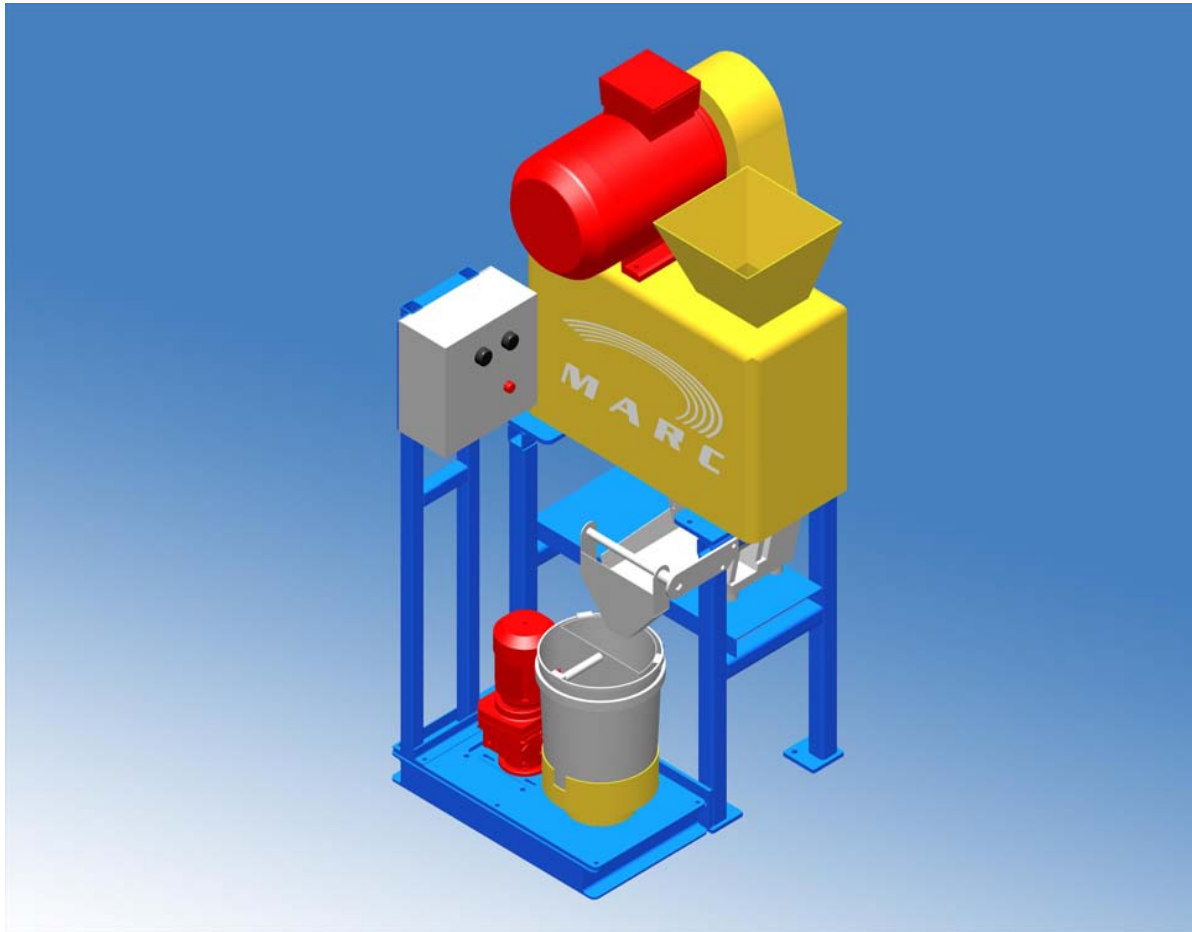

## INTRODUCE RELIABILITY INTO YOUR SAMPLE PREPARATION

MARC has combined a Jaw crusher with a bucket style RSD (Rotary Sample Divider) to produce a unique unit that is suitable for high volume sample throughput where samples are 2-3 kg or less and the residue is not required.

The rotating bucket contains a segment of your choice (generally 25%) and once the sample has been processed the segment is removed and the residue remains in the bucket.

After 10 or more samples the bucket can be removed and the collected residues disposed of into a larger bin.

Should a duplicate sample be required then additional segments can be added.

*MARC provides  
"Innovative and Specialised Engineered Products and Services"*

Tel: 61 8 9356 4000  
Fax: 61 8 9350 9333  
Email: [info@marcenv.com.au](mailto:info@marcenv.com.au)  
Web: [www.marcenv.com.au](http://www.marcenv.com.au)

